



PRODUCTS AVAILABLE
THROUGH NHS SUPPLY
CHAIN FRAMEWORK
AGREEMENT
FAG000017184

Female Health & Gynaecology

Introduction

Eurosurgical Ltd was founded in 1988 to specialise in the sales and marketing of surgical equipment, instruments and devices to the medical profession in disciplines such as Plastic and Breast Surgery, Reconstructive and Aesthetic Surgery, Gynaecology & Obstetrics, Electrosurgery and Regenerative Medicine.

Eurosurgical Ltd are the sole UK distributor for Prima Medical.

Eurosurgical Limited

Merrow Business Centre, Guildford,
Surrey GU4 7WA, United Kingdom

+44 (0) 1483 456 007
sales@eurosurgical.co.uk
www.eurosurgical.co.uk



Contents

- 03** Pipelle de cornier, Pipelle Mark II
- 04** Pipelle H
- 05** Endobrush, Scrinet
- 06** Os finder, Cervical biopsy forceps (single-use)
- 07** Cervical Dilator, Retrievette
- 08** Clear Air² Smoke (surgical plume) Evacuation system
- 09** Lletz Colposcopy loops
- 10** Electrosurgery balls and needles
- 11** WISAP Cold Coagulator and probes
- 12** WISAP C3 mobile coagulator
- 13** Cryotherapy system
- 14** Cryotherapy tips, Endosampler
- 15** Pap smear device, mucous sampler, Cytobrush
- 16** Vacuum assisted delivery cups, Amniopic, Word Bartholin catheter, Endocervical curette
- 17** Manual Vacuum Aspiration kits, Double Valve aspiration kit
- 18** Vacuum curettes
- 19** Collection sets (tubing connectors)
- 20** Cervical biopsy forceps (reusable)
- 21** HRT Tray
- 22** Cavaterm Endometrial Ablation

Pipelle

STERILE EO



Endometrial biopsy by aspiration, for effective endometrial cancer screening

- The Pipelle is the gold standard for the screening of endometrial cancers, with more than two hundred publications attesting to its effectiveness, reliability and safety.
- Its aspirating force provides high-quality sampling.
- Numbered guide-marks to gauge the insertion length.



Product name	Order Reference	NPC Code	Ø Ext.	Ø Int.	Usefull length	Ø Lateral opening	Sampling type	Packaging
Pipelle	1103000	FGR947	3.1 mm	2.6 mm	23.5 cm	2.1 mm	Histology	25 / box

Pipelle Mark II

STERILE EO



Histological and cytological sampling during the same procedure

- The Pipelle Mark II, an improved version of the Pipelle, has a rough distal tip providing for the sampling of endometrial surface cells.
- For an even more thorough sample collection, allowing for histology and cytology in just one procedure.
- Numbered guide-marks to gauge the insertion length.




← Rough tip for cytology




Product name	Order Reference	NPC Code	Ø Ext.	Ø Int.	Usefull length	Ø Lateral opening	Sampling type	Packaging
Pipelle Mark II	111020100	FGR955	3.1 mm	2.6 mm	23.5 cm	2.1 mm	Histology + cytology	25 / box


Easy to use

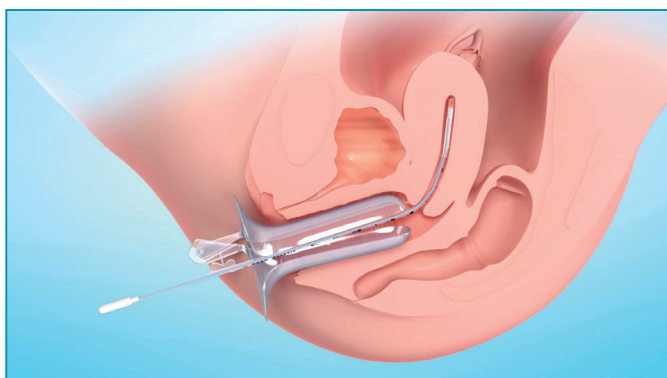

1- Insert the Pipelle until reaching the uterine fundus

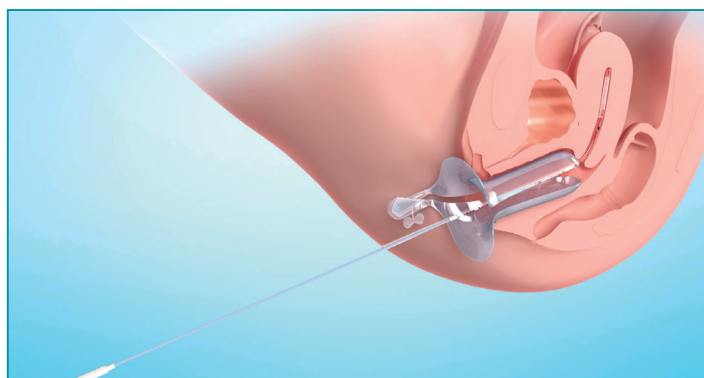

2- In one movement, pull on the plunger to bring it as far as it will go

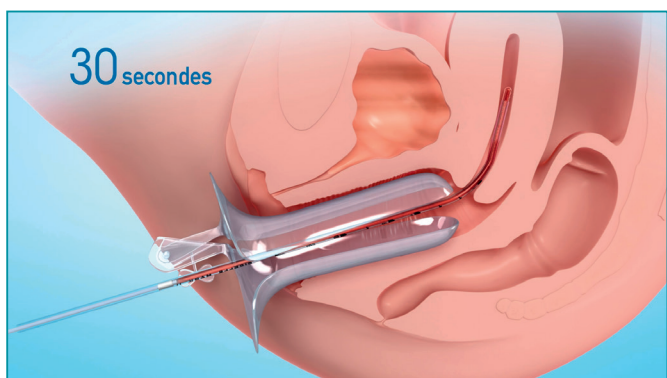

3- Sweep the uterus for 30 seconds, in back-and-forth movements and rotations



4- Cut the tip of the Pipelle



5- Push the plunger back in again to empty the content

Pipelle H STERILE EO

Based on Pipelle, Pipelle H is lengthened to allow being passed through the diagnostic sheath of a hysteroscope, once the optic has been removed.

- Owing to its length (50 cm), the Pipelle H produces enough suction to collect a large volume of endometrial tissue, despite the presence of the distension medium, and without having to repeat the procedure several times.
- As the Pipelle H is much longer than the hysteroscope, it is always possible to perform a biopsy, even if the uterine cavity is large.
- Numbered guide-marks to gauge the insertion length.
- Pipelle H is 3.1mm wide and will pass through the sheath of Hysteroscope of 3.6mm diameter



PIPELLE H 50cm 3.1mm 2.6mm 2.1mm 10/box

Product name	Order Reference	NPC Code	Ø Ext.	Ø Int.	Usefull length	Ø Lateral opening	Sampling type	Packaging
Pipelle H	1112301	FGR948	3.1 mm	2.6 mm	50 cm	2.1 mm	Histology	10 / box

Class II sterile medical devices. Manufacturer: Laboratoire CCD - CE0459. Carefully read the instructions on the labelling. Use before the expiry date indicated on the wrapping. Do not use if the wrapping or primary packaging is damaged.

Endobrush[®]


STERILE EO



Cytology and endometrial sampling

Diamond-shaped brush for sampling of endometrial surface cells. The brush retracts fully into the sheath and is covered by a rounded end cap.

- Atraumatic during passage into the cervix.
- Protection of the brush and sample against contamination.

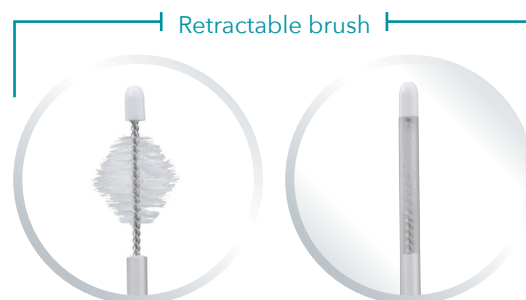

Product name	Order Reference	NPC Code	Diameter	Length	Packaging
Endobrush	1113000	FGR971	Diamond-shaped brush Length : 2 cm - Max. diameter: 15 mm	23.5 cm	25 / box

Class I sterile medical devices. Manufacturer: Laboratoire CCD - CE0459. Carefully read the instructions on the labelling. Use before the expiry date indicated on the wrapping. Do not use if the wrapping or primary packaging is damaged.


Scrinet[®]

STERILE EO



Cytology of the endocervix or bacteriological sampling

- Cylindrical polyamide swab brush, mounted on a polypropylene handle.
- Swabbing of the endocervix to sample cells for smear for cervical cancer screening or testing for bacterial infections.

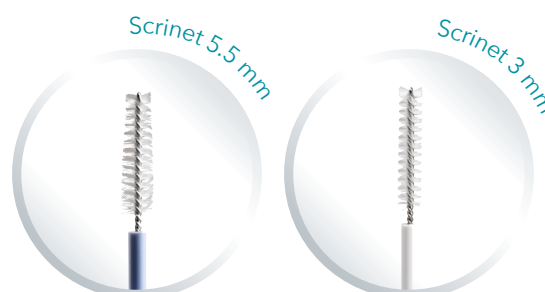

Product name	Order Reference	NPC Code	Diameter	Packaging
Scrinet 3 mm	1107030	FGR961	3 mm	50 / box
Scrinet 5.5 mm	1107055	FGR961	5.5 mm	50 / box

Available in non-sterile version (packaging: 250/box and bag of 10 units):

- Scrinet 3 mm : 1107530
- Scrinet 5.5 mm : 1107555

Scrinet: Class I sterile medical devices. Duo-Brush: Class I non-sterile medical device. Manufacturer: Laboratoire CCD - CE0459. Carefully read the instructions on the labelling. Use before the expiry date indicated on the wrapping. Do not use if the wrapping or primary packaging is damaged.

OS Finder STERILE EO

Ideal for cervical dilation, uterine sounding or other intrauterine procedures. Made of a malleable plastic and packaged individually sterile - no latex.



Product name	Order Reference	Packaging
OS Finder	022731	Pack of 50 sterile

Single Use Instruments

MedGyn is now offering disposable biopsy punches and forceps. This is the first set of instruments that are single-use.

- Tischler Biopsy Punch - Most popular bite designed for maximum tissue near os with minimal trauma
- Baby Tischler Biopsy Punch - Smaller bite size designed for less trauma and sampling
- Kevorkian Biopsy Punch - Popular design with strong narrow jaw for insertion and sampling
- Alligator Forcep - Perfect design for IUD retrieval with hinged serrated jaws



- 029701-D** Tischler Biopsy Single-Use Forceps (10/Pack)
029702-D Baby Tischler Biopsy Single-Use Forceps (10/Pack)
029703-D Kevorkian-Younger Biopsy Single-Use Forceps (10/Pack)
030112-D Alligator Single-Use Forceps (10/Pack)

Gradual Cervical Dilator

STERILE EO



No loss of cervical dilation during procedure.

The gradual dilator includes a white plastic cylindrical tube prolonged by a curved section for intra cervical insertion. The curved section is made of 4 or 5 segments (depending on the model), each being 1.5cm long and with increasing diameters as of the distal extremity. Each segment fits over the truncated end of the previous one to avoid sharp edges.

The dilators distal extremity is rounded. The proximal end includes a flat section to determine the direction of its curved section.

2 models are available; Both have a total length of 21.5cm

- Dilator Small Model (SM) dilation to 5,6,7 and 8mm
- Dilator Large Model (LM) dilation to 8,9,10,11 and 12mm



Product name	Order Reference	Packaging
Gradual Cervical Dilator SM	1605000	Box of 10 sterile
Gradual Cervical Dilator LM	1605001	Box of 10 sterile

Retrievette®

STERILE EO




IUD thread retriever

The Retrievette® is used to grab the thread or body of an intrauterine device to remove it. Polypropylene device including:

- A concave section with 3 notches and a rounded tip,
- A handle, in the lower section, to help handling.

While rotating the Retrievette® and pressing it gently against the wall of the uterus, the IUD's thread is hooked into the notches and pulled out through the cervix.

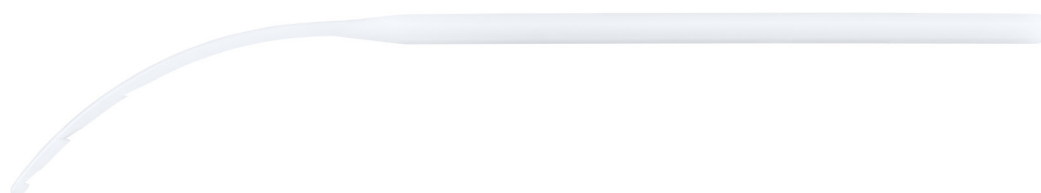

Product name	Order Reference	Packaging
Retrievette	1504000	10 / box

Prima Medical Clear Air²™

Surgical Smoke (plume) Evacuation System

Surgical plume is created by the vaporising of human tissue by surgical electro-cautery. This noxious and odorous by-product contains both organic and inorganic matter. Plume can also obscure visualisation of the tissue which could be a potential patient safety hazard. The plume falls into two categories; chemical and bacteriological, both of which have their own health risks.

The Prima Medical Clear Air² Surgical smoke evacuation system is designed to remove surgical plume particles and odour from the operating theatre, using an advanced filter system that captures 99.9% of both categories of surgical plume.

Order Reference Codes

PMS-CA1000
PMS-CAFIL


Control Unit
Replacement filters




Prima Lletz Loops and Balls




PMS1742
Single use Colposcopy
loop electrode
10mm x 10mm,
Pack of 10

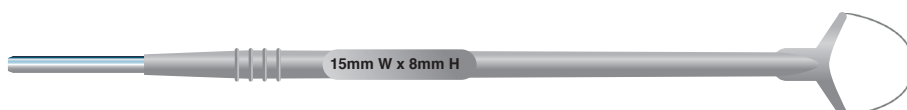

PMS1746
Single use Colposcopy
loop electrode
15mm x 8mm,
Pack of 10

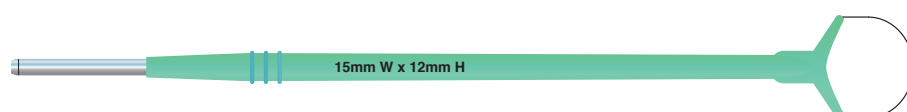

PMS1747
Single use Colposcopy
loop electrode
15mm x 12mm,
Pack of 10

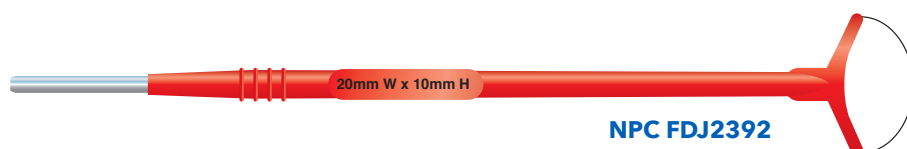

NPC FDJ2392

PMS1748
Single use Colposcopy
loop electrode
20mm x 10 mm,
Pack of 10




PMS1743
Single use Colposcopy
loop electrode
15mm x 15mm,
Pack of 10




PMS1744
Single use Colposcopy
loop electrode
20mm x 20mm,
Pack of 10



PMS1732
Single use Colposcopy
loop electrode
10mm x 8mm square,
Pack of 10

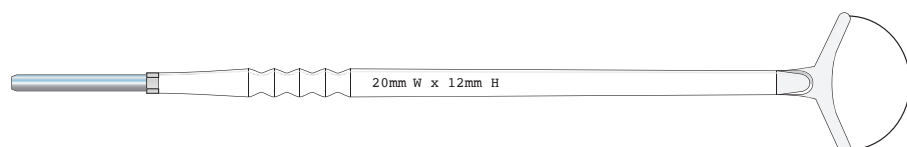

PMS1751
Single use Colposcopy
loop electrode
20mm x 12mm,
Pack of 10

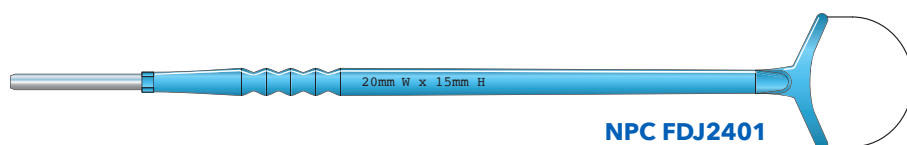

NPC FDJ2401

PMS1745
Single use Colposcopy
loop electrode
20mm x 15mm,
Pack of 10

Prima Lletz Loops and Balls are available in the standard 13.5cm length and also in a shorter length 7.0cm length to use with SuperVac Smoke Evacuation pencil with nose cone -PMS6300N

To order short Lletz Loop suffix the PMS order code with S
such as PMS1745S 20mm x 15mm loop 70mm length
PMS1744S 20mm x 20mm loop 70mm length

PMS 2905

Nonstick ball electrode
4mm ball, 50mm long, Single use,
Pack of 10

PMS 2909

Nonstick ball electrode
5mm ball, 70mm long
Single use, Pack of 10

PMS 2908

Nonstick ball electrode
3mm ball, 110mm long
Single use, Pack of 10

PMS 2907

Nonstick ball electrode
4mm ball, 110mm long
Single use, Pack of 10

PMS 2906

Nonstick ball electrode
5mm ball, 110mm long
Single use, Pack of 10

PMS 2934

Nonstick needle electrode
152mm
Single use, Pack of 10

PMS 2935

Nonstick needle electrode
152mm, shielded tip
Single use, Pack of 10

PMS 2945

Nonstick needle electrode
100mm
Single use, Pack of 10

PMS 2946

Nonstick needle electrode
100mm, shielded tip
Single use, Pack of 10

PMS 2922

Nonstick needle electrode
70mm
Single use, Pack of 10

PMS 2972

Nonstick needle electrode
70mm, shielded tip
Single use, Pack of 10

**NPC FDJ2412**

An extra long 130mm version of this ball electrode is available by adding "L" to the end of the code.



An extra long 130mm version of this ball electrode is available by adding "L" to the end of the code.

**NPC FDJ2377**

An extra long 130mm version of this ball electrode is available by adding "L" to the end of the code.

**NPC FDJ2370****NPC FDJ2367**

WISAP Cold Coagulator and Probes



APPLICATION

- HPV (CIN 1-3)
- Benign erythroplakia
- Cervix endometriosis
- Ovula Nabothi
- Chronic Cervicitis
- Hemostasis after implementation of knife conization

Cold Coagulation

Ablation of cervical lesions by „boiling“ of tissue

- Low-temperature-treatment (60-120°C)
- Fast Procedure (20 – 45 seconds)
- Treatment depth (4-7mm)
- No tissue-scarring
- Low pain scores
- Low bleeding events
- Low watery discharge

Mains connection:	Standard: 230V - 50 /60 Hz Please indicate other requested voltages when ordering
Power Input:	66VA
Range of temperature:	0-120 ° C
Regulation of temperature:	Stepless
Dimensions:	W x H x D = 270 x 120 x 190 mm
Casing:	Stainless steel, aluminium/white coated
Weight:	3,6 kg
Order data	CERVIX COAGULATOR
6001	Cervix-Coagulator
6020	Instrument cable, orange, with two coupling plugs, silicone, L = 1,5 m
600 - X	Therapy Probes


TREATMENT DEPTH

at 80 -100°C

5" 1mm

10" 2mm

15" 3mm

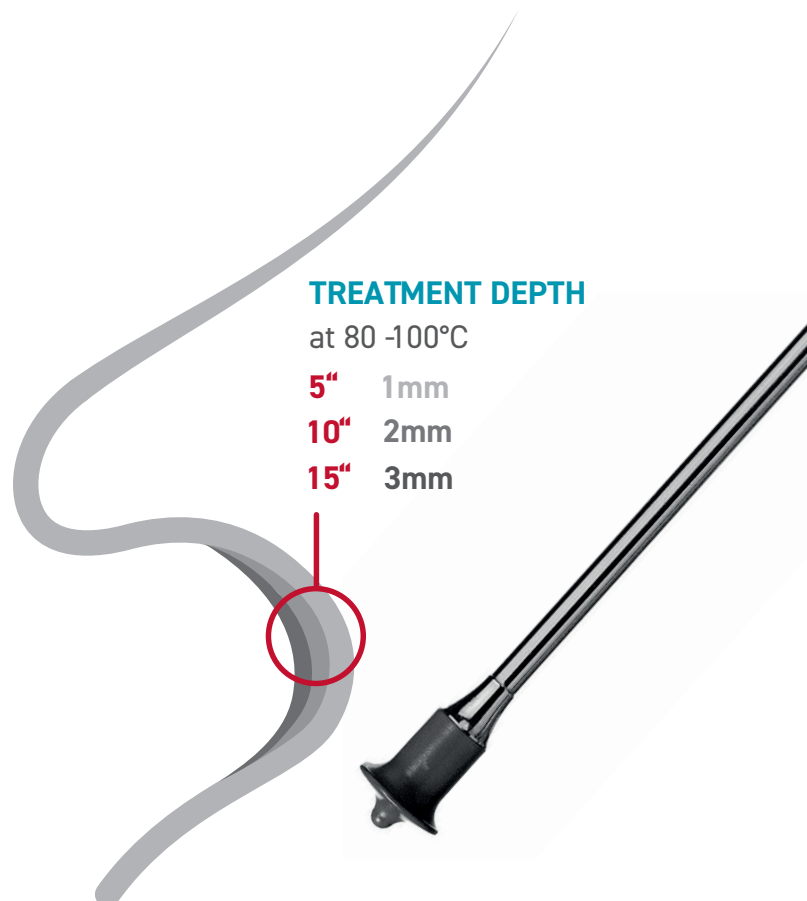



Model A Ref 6002; For large area erythro and leucoplakia



Model E Ref 6006; For endo-uterine therapy 120mm



Model B Ref 6003; For ectopies at dehiscant cervical os



Model F Ref 6007; For large areas of erythron and leucoplakia



Model C Ref 6004; For small circular ectopies



Model G Ref 6008; For ectopies at dehiscant cervical os



Model D Ref 6005; For endo-cervical therapy, max 60-70°C



Model H Ref 6009; For wounds of conisation

WISAP C3 Mobile Coagulator

SIMPLE HANDLING for effective treatment of cervical precancer

Globally, there is a pressing need for a simple treatment for cervical precancer. Current methods are expensive and difficult to implement on a large scale: This LMIC-adapted thermo coagulator is designed with one idea in mind: Simple handling for maximum effect. Therefore the device features a single optimally designed thermo-probe, together with a fixed temperature setting and a timer function. That way we eliminate any possible source for human error and ensure the greatest efficacy on site.



Cryotherapy System

MedGyn Cryotherapy System is designed by physicians and engineered by cryogenic experts to ensure safety, performance and convenience. The MedGyn Cryotherapy System has superior technology combined with a 3-position single trigger to deliver optimal cryosurgical temperatures. A wide variety of autoclavable tips is available. The equipment works with Nitrous Oxide (N₂O) or with Carbon Dioxide (CO₂) gas and has standard- specific adaptors.



019000 MedGyn Cryotherapy System MGC-200 (tips not included)

- Single-hand control from three-position trigger (freeze, off, defrost)
- Instant defrost
- Wide choice of autoclavable tips
- Valve body designed and manufactured for long, trouble-free life
- Ability to change tip during procedure without shutting off gas flow
- Capture "O" ring design provides complete gas seals where tips attach to probe system
- Built-in regulators control pressure at tip for added safety and gas economy
- Nitrous Oxide (N₂O) and Carbon Dioxide (CO₂) adaptors available

Cryotherapy Equipment & Accessories

- 019019 Carrying Case for Cryotherapy System
- 019020 Protective Cap for Cryotherapy System
- 019021 Rubber Plug 6448K92 for MedGyn Cryo Tip
- 019030 450 Purifier Housing with 1/4 npt Male Thread Adapter 019031 541 Filter Cartridge for Purifier
- 019011 British Adapter for N₂O
- 019001 Cart for Cylinder, 6lb or 20lb 019003 «O» Ring for Cryo System (3/Pack) 019008 French Adapter for CO₂

Cryotherapy Tips



019002/2507 MedGyn Cryotherapy Tip GT-2507 Endo/Exocervical Large



019002/1910 MedGyn Cryotherapy Tip GT-1910 Endo/Exocervical Small



019002/1900 MedGyn Cryotherapy Tip GT-1900 Exocervical, 19mm Flat



019002/2500 MedGyn Cryotherapy Tip GT-2500 Exocervical, 25mm Flat



019002/1905 MedGyn Cryotherapy Tip GT-1905 Exocervical Convex (19mm x 5mm)



019002/0811 MedGyn Cryotherapy Tip GT-0811 Skin Lesion, 8mm Diameter 45° angle



019002/0800 MedGyn Cryotherapy Tip GT-0800 Endocervical, Round



019002/0519 MedGyn Cryotherapy Tip GT-0519 Endocervical (Nulliparous)



019002/0507 MedGyn Cryotherapy Tip GT-0507 Skin Lesion, 5mm Diameter 45° angle



019002/0219 MedGyn Cryotherapy Tip GT-0219 Micro, 2mm Diameter

MedGyn Endosampler

Utilized for sampling of endometrial tissue, the unique design of the MedGyn Endosampler™ and Endosampler II™ allows for superior accuracy and a more reliable sample than traditional endometrial sample devices.

The Endosampler II™ is available with two holes for easier and more comfortable suction.

- Includes 3mm curette & 10cc syringe with lock-spring mechanism
- Highly-visible, graduated markings on curette
- Semi-rigid curette allows for easier insertion and injection of saline or contrast
- Curved design for easier manipulation
- Back-flow prevented by syringe lock-spring mechanism
- Sharper slot at curette tip for better sampling
- "Closed" system prevents sample contamination

022720 MedGyn Endosampler™ (25/Box)

022722 MedGyn Endosampler II™ (25/Box)

NPC Code FGR949



MedGyn Mucous Sampler

The MedGyn Mucous Sampler is a single-use, sterile, endocervical aspirator.



- The device consists of a 23 cm long, clear, flexible polypropylene tube.
- The distal tip of the device is rounded with a circular opening of 2.2 mm diameter for the specimen to be aspirated into the sheath lumen.

022729 Disposable Mucous Sampler single-use Pack of 25

MedGyn Pap Smear Device

The MedGyn Cell Sweep® pap smear device simultaneously samples both the endocervix and ectocervix eliminating the need for both a Cytobrush and a Spatula.



- Allows for optimal cell harvesting and cell transferring
- Dual sampling of endo and ectocervix and separation of cells in one-slide technique
- Minimizes false negatives
- Minimizes red blood cells and blood artifacts

022362 MedGyn Cell Sweep® Slide Holder (100/Pack)

022361 MedGyn Cell Sweep® Non-Sterile (100/Pack)

022361K MedGyn Cell Sweep® Non-Sterile with Kit (Slide and Slide Holder) (100/Pack)

022363 MedGyn Cell Sweep® Sterile (100/Pack)

022363K MedGyn Cell Sweep® Sterile with Kit (Slide and Slide Holder) (100/Pack)

MedGyn Cytobrush


Standard devices for one of the most frequently performed gynecology procedures.



022710 MedGyn Cytobrush (100/Box)
022364 MedGyn Pap Cell (50/Box)

Medgyn vacuum assisted delivery cups

MedGyn Hand-Held Vacuum Delivery Pump may be used to safely assist medical personnel during vaginal deliveries. MedGyn Vacuum Delivery Cups offer optimal adhesion to ensure proper vacuum levels.

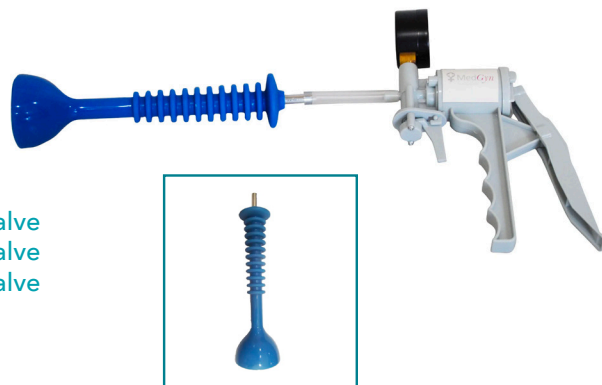

- 020322** MedGyn Vacuum Delivery Cups 50mm with Release Valve
- 020323** MedGyn Vacuum Delivery Cups 60mm with Release Valve
- 020316** MedGyn Vacuum Delivery Cups 70mm with Release Valve
- 020320** Vacuum Delivery Cup 50mm
- 020321** Vacuum Delivery Cup 60mm **NPC Code FGR969**
- 016160** MedGyn HVP-100 Hand-Held Vacuum Delivery Pump

Amniopic

Ergonomic and safe amniotic membrane perforator.

- Slightly curved amniotic membrane perforator.
- Rounded exterior and profiled hook to prevent lesion of surrounding tissues during insertion.



16004PO Amniopic 25 cm (25/pack) **NPC Code FGR966**



MedGyn Word Bartholin Catheter

For use following an incision and drainage of a Bartholin's cyst or abscess. Catheter kit includes balloon catheter, syringe for balloon inflation and scalpel.



022719 MedGyn Word Bartholin's Catheter (6/Box)

NPC Code FPF278

MedGyn Endocervical Curette

The MedGyn Endocervical Curette provides an efficient sampling of the cervical canal to test for abnormal, precancerous conditions, or cervical cancer.

- The Kevorkian-Young endocervical curette is a sterile, single-use device with a basket and a square tip.

022732 MedGyn Disposable Endocervical Curette (25/Box)



Medgyn MVA Aspiration kits

MedGyn Aspiration Kit and Aspiration Kit IV are effective devices for simultaneous endometrial curettage and aspiration.

- 022514** MedGyn Aspiration Kit with 4mm Cannula-Curette (10/Pack) **NPC Code FPF282**
- 022515** MedGyn Aspiration Kit with 5mm Cannula-Curette (10/Pack) **NPC Code FPF283**
- 022516** MedGyn Aspiration Kit with 6mm Cannula-Curette (10/Pack) **NPC Code FPF280**
- 022517** MedGyn Aspiration Kit with 7mm Cannula-Curette (10/Pack) **NPC Code FPF281**
- 022511** MedGyn Aspiration Kit IV with 4mm, 5mm, 6mm and 7mm Cannula-Curette (10/Pack)



MedGyn Double Valve Aspiration

MedGyn Aspiration Kit Double Valve is used for uterine aspiration/ uterine evacuation in obstetrics and gynecologic procedures.

The device fits the MedGyn Flexible Karman Curettes and a wide variety of MedGyn Rigid and Straight Curettes.



- 022526** Sterile, single-use packaged individually **NPC Code FPF279**

Medgyn Vacuum Currettes

MedGyn offers a complete line of high-quality, disposable vacuum cannula-curettes in the industry. Available in Rigid Curved, Rigid Straight, Flexible and Flexible Karman options, our cannula-curettes are individually packaged and sterilized.

Vacuum currettes are available with standard ½" (12.7mm) or ⅜" (9.52mm) hubs, that are designed to fit the appropriate Collection Set (see page 19). Please note, Collection Set 022310-WO DOES contain an adaptor that converts from ½" to ⅜" to accommodate both hub sizes of currettes.

All available sterile, single-use packaged individually

Rigid Curved Cannula-Curette

022106	6mm NPC Code FGR956
022107	7mm NPC Code FGR952
022108	8mm NPC Code FGR962
022109	9mm
022110	10mm NPC Code FGR951
022111	11mm NPC Code FGR964
022112	12mm NPC Code FGR967
022113	13mm - ⅜" Base
022114	14mm - ½" Base
022144	14mm - ½" Base w/Adapter
022145	14mm - ⅜" Base
022115	15mm - ½" Base
022116	16mm - ½" Base



Rigid Straight Cannula-Curette

022206	6mm
022207	7mm
022208	8mm
022209	9mm
022210	10mm NPC Code FGR958
022211	11mm NPC Code FGR964
022212	12mm NPC Code FGR954



Flexible/Karman Cannula-Curette

022213	13mm - ⅜" Base
022214	14mm - ½" Base
022244	14mm - ½" Base w/Adapter
022245	14mm - ⅜" Base
022215	15mm - ½" Base
022216	16mm - ½" Base

Flexible Cannula-Curette with Dual Port Design

022004	4mm
022005	5mm
022006	6mm
022007	7mm NPC Code FGR970
022008	8mm
022009	9mm
022010	10mm
022011	11mm
022012	12mm



Flexible/Karman Cannula-Curette

022004K	4mm
022005K	5mm
022006K	6mm
022007K	7mm
022008K	8mm
022009K	9mm
022010K	10mm
022011K	11mm



Medgyn Collection Sets

MedGyn Disposable Collection Sets include 6' tubing (either 3/8" or 1/2" I.D.) attached to a control handle with slip-ring for finger-tip vacuum control. The Collection Sets are individually packaged and sterile.

Collection Sets & Bottles

022310	3/8" Collection Set with Handle
022310NS	3/8" Collection Set w/Handle & Adapter (Non-Sterile)
022310-WO	3/8" Collection Set with Handle (Without Adapter, Sterile) (50/Pack) NPC Code FGR968
022310-50	3/8" Collection Set with Handle (50/Pack) NPC code FGR953
022314	3/8" Collection Set x 8 Ft with Handle & Adaptor
022312	1/2" Collection Set with Handle & Adapter
022312-25	1/2" Collection Set with Handle (25/Pack)
022312NS	1/2" Collection Set w/Handle & Adapter (Non-Sterile)
022331	3/8"Disp. Collection Jar - with Cap (10/Box)
022332	1/2" Disp. Collection Jar - with Cap (10/Box)
003984-901	3/8" Disposable Collection Jar with Cap & Tissue Trap (10/Box)
022323	Disp. Tops Large Non-Conductive (10/Box)
022324	Disp. Tops Small Non-Conductive (10/Box)

Accessories

022316-3/8-6	Tubing 3/8" Precut to 6'
022311	Tubing 3/8" Custom Cut Per Foot
022316-1/2-6	Tubing 1/2" Precut to 6'
022313	Tubing 1/2" Custom Cut Per Foot
022317	Tubing 16" Connector
022318	Tubing 20" Connector
022317-ADAP	Tube Fitting - Male 3/8" Plastic
022317-INT	Tube - Interconnecting with Adapter (10/Box)
022350	Plastic Adapter for 1/2" Base Curettes (Sterile)
022351	Plastic Adapter for 1/2" Base Curettes (Non-Sterile)
O-RING-F.AD	3/8" O Ring for Female Adapter (10/Pack)
022320	Disposable Filter Assembly (10/Box)
T.TRAP/ST	Safe-Touch Tissue Trap (10/Box)
ADAP-REDC	1/2" to 3/8" Adapter for Aspiration Handle (For Small Curette) (10/Box)
022354	3/8" to 1/2" Adapter for Aspiration Handle (For Large Curette)



Permeable Gauze Sacks

Gauze collection sack easily attaches to collection canister for simple and efficient separation of tissue from aspirated fluid.


022460	Gauze Sacks Tissue Trap (10/Box)
O-RING3/8	O Ring for 3/8" Sack (10/Box)
O-RING1/2	O Ring for 1/2" Sack (10/Pack)

Metal Swivel Handles - Stainless steel, autoclavable, with slip-ring control

022450	3/8" Handle for Curettes up to 14mm with small base
022451	1/2" Handle for 14mm, 15mm & 16mm Curettes
022340	Stainless Steel 3/8" Female Fitting
022341	Stainless Steel 3/8" Male Fitting
022343	Stainless Steel 1/2" Male Fitting
022352	Stainless Steel 3/8" Handle Extender
022353	Stainless Steel 1/2" Handle Extender

Reusable Instruments for Cervical biopsy

Rotating Cervical Biopsy Punch Forceps

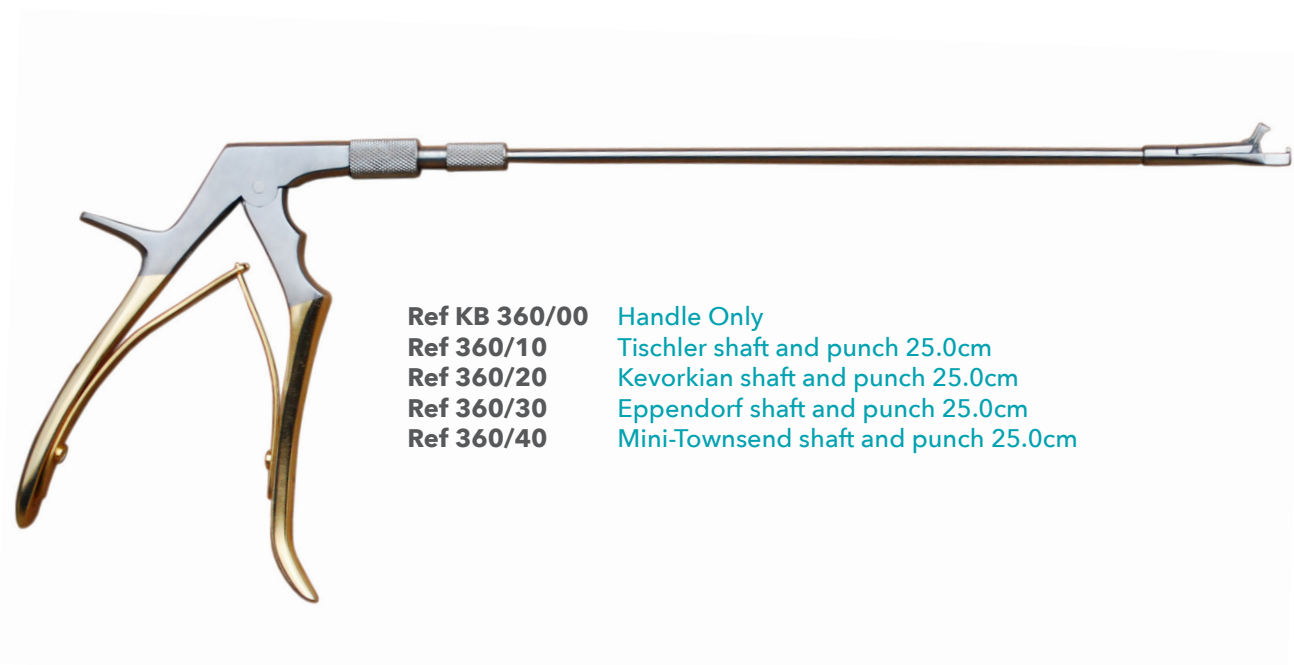

Non-Rotating Cervical Biopsy Punch Forceps – handle and shaft/punch combined



A wide range of other Gynaecology reusable instruments are available on request, including;

Double-Ended Dilators (Pratt, Hegar and Hank styles)
 Sims Curettes
 Endocervical Curettes
 Uterine sounds
 Ovum forceps
 Myoma Grasping forceps
 Allis Tissue Forceps
 Placenta Forceps
 Polyp Forceps
 Hysterectomy Clamps
 Metal Needle Extenders (luer lock hubs)

HRT Tray

A sterile single-use tray / kit designed for placement of Hormone Replacement Therapy (HRT) is available, to help nurses and doctors to reduce time when performing hormone implant therapy treatments. The HRT kits eliminates both expensive cleaning, decontamination and sterilisation costs and provide clean, sterile and new instruments for every patient.

The HRT kit Contains all the necessary equipment for hormone implantation except the implants.



This HRT Pack consists of the following:

- 1 x trocar and cannula
- 1 x B/P handle with 15 blade
- 1 x 5ml syringe luer slip
- 1 x green needle 21 G x 1.5 luer slip
- 1 x blue needle 23 G x1.25 luer slip
- 1 x Addson dissecting non toothed forcep
- 5 x non woven Gauze swabs 10x10cm
- 1 x fenestrated drape 76x76 cm with window 10x10cm
- 1 x tray 2 compartments Poly with integrated gallipot
- 1 x trace me lable
- 5 x Balls non woven size Medium

Order Reference Code GPSS327

Cavaterm Endometrial Ablation

Cavaterm™ is a minimally invasive treatment for Dysfunctional Uterine Bleeding (DUB)

A safer and more flexible endometrial ablation system, that is effective in over 90% of uterine cavities with extremely high satisfaction rates when compared to alternative treatments.

The new second generation Cavaterm endometrial ablation technique consists of an adjustable, flexible silicone balloon catheter and a Cavaterm control unit.

The total treatment combines heat and pressure during 10 minutes to irreversibly coagulate the endometrium and underlying myometrium to a depth of 5-8 mm. Cavaterm treatment can be carried out in outpatient clinic and as a day case procedure.



Specific features are built in the system to ensure safety, efficacy and comfort.

The Cavaterm™ system has been successfully used since 1993. It is proven to be safe, efficient and easy to use. Over the past years further product developments based on users feedback were implemented for an optimal evolution of the system. Many thousands of patients have been treated successfully with the Cavaterm™ system. Clinical studies report average success and satisfaction rates exceeding 90%, which remain stable throughout the years.^{1,2,3}

TECHNICAL DATA

Catheter

Shaft/cervical diameter	6 mm
Balloon adjustment	4-8 cm
Target temperature	78°C
Operating pressure	230-240 mmHg

Power supply

Mains connection or rechargeable	100-240 V: 50-60 Hz
Lead batteries	24 V / 2 Ah
Full recharge time	Typically 4 hours
Battery top up	Approx. 45 minutes

Central unit

W x L x H	450 x 310 x 120 mm
Weight	8 kg

Order Reference Code for Cavaterm Catheter;
Order Reference Code for Cavaterm Control Unit;

CAV2010-10
CAV2020

References

1. Hawe J.A., Phillips G., Chien P., Erian J. & Garry R. Cavaterm thermal balloon ablation for the treatment of menorrhagia. Br J Obstet Gynaecol 1999;Nov;106;11;1143-1148
2. Friberg B. & Ahlgren M. Thermal balloon endometrial destruction: the outcome of treatment of 117 women followed up for a maximum period of 4 years. Gynaecological Endoscopy 2000; 9:389-395
3. Mettler L. Long-term results in the treatment of menorrhagia and hypermenorrhea with a thermal balloon endometrial ablation technique. JSLS 2002 6;4;305-9

Available through NHS SC Framework Agreement Minimally Invasive Surgery 2020/S 170-412525



Eurosurgical Limited

Merrow Business Centre, Guildford,
Surrey GU4 7WA, United Kingdom

+44 (0) 1483 456 007
sales@eurosurgical.co.uk
www.eurosurgical.co.uk